

535706

Nucleotide

PubMedNucleotideProteinGenomeStructurePopSetTaxonomyOMIMBooks

Search  for

LimitsPreview/IndexHistoryClipboardDetails

☐ 1: L27092. Pneumocystis cari...[gi:535706]


Related Sequences, Protein, PubMed, Taxonomy

LOCUS PMCMST1 3363 bp DNA linear PLN 26-SEP-1994  
DEFINITION Pneumocystis carinii B-cell receptor (msgI) gene, 3' end.  
ACCESSION L27092  
VERSION L27092.1 GI:535706  
KEYWORDS B-cell receptor.  
SOURCE Pneumocystis carinii DNA.  
ORGANISM Pneumocystis carinii  
Eukaryota; Fungi; Ascomycota; Pneumocystidomycetes;  
Pneumocystidaceae; Pneumocystis.  
REFERENCE 1 (bases 1 to 3363)  
AUTHORS Garbe,T.R. and Stringer,J.R.  
TITLE Molecular characterization of clustered variants of genes encoding  
major surface antigens of human Pneumocystis carinii  
JOURNAL Infect. Immun. 62 (8), 3092-3101 (1994)  
MEDLINE 94314421  
PUBMED 7518806  
FEATURES Location/Qualifiers  
source 1..3363  
/organism="Pneumocystis carinii"  
/db\_xref="taxon:4754"  
gene 152..3244  
/gene="msgI"  
CDS <152..3244  
/gene="msgI"  
/codon\_start=1  
/product="B-cell receptor"  
/protein\_id="AA21645.1"  
/db\_xref="GI:535707"  
/translation="VARAVKRQVAGVKNNEAEERLFALIMRADYKDESKCKNKIKEYC  
DGLKNASLTSEEVHKELKDFCKDGSQGGKCEELKKNVEAKCNFPKTKLEGLVKKDASG  
LTNDDCKENERQCLFLEGACPDLVEDCSKLRNLQYQKKREGVAEEVLLRALRGDLGNK  
TECEKKIKDVCPKIGQESDELTLCLDQKKCTNLMTARDKKCNTLEEDVKALENKN  
NLLGKCLPLLEHATFTEGTAKKASQCTPNKDCEDYLPKCELAEECGKGIYIHPGP  
DFDPTKPEPTVAEDIGLEELYKKAEDGVHIGKPPVRDATAALLLIQNPDPKIQANE  
KEKCKVLENKCKELKKHEVLGDLNQNAASQSGTKKCEELEKELANSTKILSEKIKN  
KHLSGSGETIPWYKLSTFLSDSDCARLESDFYFAQDKDPLKKECKNVKAACYKRGLD  
ARANKVLQENMRGLLRGSNQSWLKKFQQELVKVCEKLKEENKGSFSNDELFLVLCVQPA  
KAARLLTHDLRMTTIFLRQQLDQKRDFTVKTAKGIREKCQDLGKGFGKEITWPCHTL  
EQQCNRGLGTTIILKQVLLNEHKDTLKTENCVTYLKEKCNKWSRRGDDRFVFCVFQ  
ATCKLMVKDVQDRCKIFKENIKVSEIVDFLKNNTNNITTLERNCPSWHTYCNRFSSNC  
PDFSKNPCTKIKNNCKPFYERKALEDALKVELRGKLSDENKCTAALKGYCTLAGNVN  
NASVRSCLKDNTQGSNNKTDEKVVEELCKKLMEEVKEQCETLPAELKQPADDLEKDV  
TYEELKEEAKKAMNKSSLVLSFVKKDGNNTPKNNSKSEDKNVVSNEKDTIKHVILRR  
GVKDVLVTELEAKAFDLAAEVFGRYVDLKERCEKLTLDGKIDCDGLKGVCGKIKKK  
CRDLKPLEVKSHEIVTESTTTTTTTTTTTTVDPKATECKSLQTTDTWVTQTSTHTSTST  
ITSTITSKITLTSTRCKPTKCTTGDDAEDVKPSEGLRVSGWNVMRGAIVAMVISFMI  
"  
mat\_peptide 152..3241  
/gene="msgI"  
/product="B-cell receptor"  
/standard\_name="major surface antigen"  
/evidence=experimental

BASE COUNT	1312 a	505 c	753 g	793 t
ORIGIN				
1	gtattgtggt	gtaaaaacta	ccttttaaac	tttcaatatt
61	tataaaatag	aaacagcttt	aaaagtctgc	ccacacaccc
121	cattaatgtc	ctctagat	tggcgcgggc	ggtggcgcgg
181	aggggtaaag	aataatgaag	ctgaagaacg	tctttttgct
241	taaagatgaa	agtaaatgta	aaaataaaat	aaaagaatat
301	ttcactaaca	tcagaagaag	tacataaaaga	ttttgcaaag
361	aggaaaaaaa	tgcaagaac	taaaaaaa	tggtgaagca
421	aaaattggaa	ggattggtga	aaaaagatgc	ctcaggcttg
481	gaatgaacga	caatgcctat	ttttggaggg	agcatgccct
541	caaactaagg	aatctatggt	accagaaaaa	acgtgaagga
601	gagggcactt	cgtggggatt	taggaaataa	aactgaatgc
661	ttgccccaaa	ataggccaag	aaagtgatga	ggtgacattg
721	aacatgcacg	aatcttatga	ctgcgcgaga	caaaaagtgt
781	taaaaagca	cttgaaaaca	aaaacaactt	gttaggaaaa
841	tgctactttc	acagagggga	ctgcaaaaaa	agcatcacag
901	cgaggactat	ctaccaaagt	gtgatgaatt	ggcagaagaa
961	ttacatacat	ccgggacccg	attttgatcc	aactaagcca
1021	catagggctg	gaagagcttt	ataaaaaggc	cgcagaagat
1081	tcctgtaaga	gacgcaactg	ctctactggc	acttttgatt
1141	tcaagcgaac	gaaaaagaaa	aatgcaaaaa	agttcttgaa
1201	aaaacatgag	gttttgggag	atctttgtaa	tcaaaatgca
1261	aaagtgtgaa	gaactagaga	aggaattagc	aaacagtact
1321	aaagaataaa	catctctctg	gatccggaga	aaccattcca
1381	tcttagtgac	agtgactgtg	caaggttaga	atcagactgt
1441	agatcctctt	aaaaaagaat	gcaagaatgt	taaggcagca
1501	tgacggggca	aataaagtgc	tgcaagaaaa	tatgctggg
1561	aagttggctt	aagaaatttc	aacaagaatt	agtaaaagta
1621	aaataaagga	agtttttcaa	acgatgaatt	atgtgttctg
1681	agcccgggtg	cttacacatg	atcttcgaat	gacaactatc
1741	tcaaaagcga	gattttccga	cagtaaagac	tgccaaggga
1801	tattgggaaa	ggattccaaa	aagaaattac	gtggccatgt
1861	caatcgcttg	gggactacag	aaattttgaa	gcagggttta
1921	tttgaaaact	catgaaaact	gtgtaacata	tttaaaggaa
1981	aaggggcgat	gaccgtttct	cttttgatg	tgttttccaa
2041	ggtaaaagac	gtacaagata	ggtgcaaaat	attcaaagaa
2101	tgttgatatt	cttaaaaaata	atacgaataa	cataacaaca
2161	ctggcatacg	tactgcaata	gattttcatc	taattgtcca
2221	ttgtacaaaa	attaaaaaca	attgtaagcc	attttatgaa
2281	tctcaaagta	gagcttcgag	gaaagctgag	tgatgaaaac
2341	tgagatattg	acactagcgg	gaaatgtaaa	taacgcatac
2401	taacactcaa	ggtagcaata	aaaagactga	cgaaaaagtg
2461	attgatggaa	gaagtaaaag	agcaatgcga	gacattacca
2521	agacgatcta	gaaaaagatg	ttaagacata	tgaggaactt
2581	aatgaacaag	tccagccttg	ttttatcatt	cgtttaagaa
2641	aaataatagc	aaaagcgaag	ataagaatgt	cgtttcaaat
2701	tgtgaaaata	ctacggagag	gagttaaagg	tgtattagta
2761	gtttgactta	gcagcagaag	tgtttggcag	atatgtagat
2821	attgaccttg	gattgtggga	ttaaagacga	ttgcatgggt
2881	gattaagaag	aaatgtcgcg	atctgaagcc	tctggagggt
2941	agaaagcaca	acgacgacca	caacgacaac	aacgaccgtt
3001	atgcaaatcc	ttacagacaa	cagacacatg	ggttacacag
3061	atctaccatc	acatctacga	ttacatcaaa	aataacattg
3121	accaaccaag	tgtacgacag	gggatgatgc	agaagacgtg
3181	ggtgagcggg	tggaatgtga	tgaggggggc	aatagtagca
3241	ttagaaagca	aatgagcaag	atttaattgat	gagcatgtat
3301	ttaggacagg	aaagagtaaa	gtagattgag	tcattctctc
3361	aca			

//

Revised: October 24, 2001.

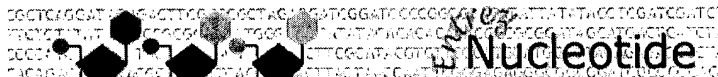

PubMed	Nucleotide	Protein	Genome	Structure	PopSet	Taxonomy	OMIM	Books
Search		Nucleotide	for				Go	Clear
Limits		Preview/Index		History		Clipboard		Details
Display	default	Save	Text	Add to Clipboard				

☐ 1: U34351. *Pneumocystis carinii* [gi:1144496]

Protein, PubMed, Taxonomy

LOCUS PCU34351 306 bp DNA linear PLN 25-APR-1996  
 DEFINITION *Pneumocystis carinii* human isolate glycoprotein A gene, partial cds.  
 ACCESSION U34351  
 VERSION U34351.1 GI:1144496  
 KEYWORDS .  
 SOURCE *Pneumocystis carinii*.  
 ORGANISM *Pneumocystis carinii*  
 Eukaryota; Fungi; Ascomycota; Pneumocystidomycetes;  
 Pneumocystidaceae; *Pneumocystis*.  
 REFERENCE 1 (bases 1 to 306)  
 AUTHORS Wright, T.W., Gigliotti, F., Haidaris, C.G. and Simpson-Haidaris, P.J.  
 TITLE Cloning and characterization of a conserved region of human and rhesus macaque *Pneumocystis carinii* gpA  
 JOURNAL Gene 167 (1-2), 185-189 (1995)  
 MEDLINE 96144272  
 PUBMED 8566775  
 REFERENCE 2 (bases 1 to 306)  
 AUTHORS Wright, T.W.  
 TITLE Direct Submission  
 JOURNAL Submitted (21-AUG-1995) Terry W. Wright, Microbiology and Immunology, University of Rochester School of Medicine and Dentistry, 601 Elmwood Avenue, Rochester, NY 14642, USA  
 FEATURES  
 source Location/Qualifiers  
 1..306  
 /organism="Pneumocystis carinii"  
 /specific\_host="Homo sapiens"  
 /db\_xref="taxon:4754"  
 /clone="Pc Gibbs 300"  
 /clone\_lib="PCR Generated Fragment"  
 CDS  
 <1..>306  
 /note="gpA"  
 /codon\_start=1  
 /product="glycoprotein A"  
 /protein\_id="AAA97603.1"  
 /db\_xref="GI:1144497"  
 /translation="RVACYKKGLDAQANEVLQSKMHELFRGSGEEWFKKLLDKITKEC  
 QKLKTASDELFLLCIEPFIIVRILAADFRLRMVFLQEQNLNENRDFPTKENCVELKEKC  
 "  
 primer\_bind 1..23  
 /note="Cys-Sense primer"  
 primer\_bind complement(283..306)  
 /note="Cys-Antisense primer"  
 BASE COUNT 116 a 43 c 68 g 79 t  
 ORIGIN  
 1 agagtagcgt gttataagaa gggacttgac gcacaggcaa acgaagtatt acagagtaaa  
 61 atgcacgaat tggtccgtgg ttctggagaa gaatggttta aaaaattact ggataaaata  
 121 acaaaagaat gtcaaaagct caaaacagca agtgacgagt tattttctact atgcattgaa  
 181 ccatttatag cggttcgaat acttgcagca gattttcgat taagaatggt ttctctacaa  
 241 gaacaactaa atgaaaaccg ggattttcct acaaaagaaa actgtgttga gttgaaggag  
 301 aagtgt  
 //

NCBI PubMed Pub? QUERY PubMed ?

Other Formats: Citation MEDLINE

Links: Related Articles

☐ Order this document*Mol Cell Probes* 1991 Apr;5(2):143-149

## The use of europium (Eu<sup>3+</sup>) labelled primers in PCR amplification of specific target DNA.

Dahlen P, Iitia A, Mikkala VM, Hurskainen P, Kwiatkowski M

Pharmacia Genetic Engineering Inc., La Jolla, CA 92037.

The polymerase chain reaction (PCR) has many potential applications in the field of DNA probe diagnostics. Here we describe a method that utilizes PCR and time-resolved fluorometry (TRF) for the detection of specific target DNA. First the DNA segment to be detected is amplified according to standard procedures. Then a pair of europium (Eu<sup>3+</sup>) and biotin-labelled primers nested within the amplified fragment is incorporated in a few additional PCR cycles. Thus amplified DNA fragments are generated that contain an affinity label (biotin) and a detectable label (europium). The doubly-labelled amplified DNA fragments are collected onto streptavidin coated microtitration strips and the bound Eu<sup>3+</sup> is measured in a time-resolved fluorometer. We show here the application of this method to the detection of HIV-1 DNA. As few as five copies of HIV-1 DNA could readily be detected using this assay. The method described here is sensitive, rapid and easy to employ. In addition it lends itself to automation.

PMID: 2072935, UI: 91304446

Save the above report in Macintosh Text format

Order documents on this page through Loansome Doc

NCBI PubMed Pub MEDLINE QUERY PubMed ?

Other Formats: Citation MEDLINE

Links: Related Articles

☐ Order this document*J Clin Microbiol* 1991 Apr;29(4):798-804

## Detection of human immunodeficiency virus type 1 by using the polymerase chain reaction and a time-resolved fluorescence-based hybridization assay.

Dahlen PO, Iitia AJ, Skagius G, Frostell A, Nunn MF, Kwiatkowski M

Pharmacia Genetic Engineering Inc., La Jolla, California 92037.

The polymerase chain reaction (PCR) has many potential applications in the field of nucleic acid diagnostics. In particular, it has been successfully applied to the detection of pathogens present in low copy numbers such as the human immunodeficiency virus type 1. Here we describe a time-resolved fluorescence-based hybridization assay which, combined with the PCR, offers an extremely sensitive method for the detection of nucleic acids. In this assay format, the PCR is run by standard procedures and the subsequent hybridization reaction is carried out in solution by using two oligonucleotide probes, one biotinylated and one labeled with europium (Eu3+). The sandwich hybrids are then collected onto a streptavidin-coated microtitration well, and the bound Eu3+ is measured in a time-resolved fluorometer. This assay is rapid, user friendly, and quantitative and lends itself to automation. The application of this assay to the detection of human immunodeficiency virus type 1 is described.

PMID: 1890180, UI: 91365881

 the above report in   format documents on this page through Loansome Doc

NCBI PubMed PubM QUERY PubMed ?

Other Formats: Citation MEDLINE

Links: Related Articles

☐ Order this document*Anal Chem* 1993 Jan 15;65(2):158-163

## Quantification of polymerase chain reaction products in agarose gels with a fluorescent europium chelate as label and time-resolved fluorescence spectroscopy.

Chan A, Diamandis EP, Krajden M

Department of Clinical Biochemistry, Toronto Hospital, Ontario, Canada.

We have 5'-end-labeled one polymerase chain reaction (PCR) primer with the europium chelator 4,7-bis(chlorosulfonylphenyl)-1,10-phenanthroline-2,9-dicarboxylic acid (BCPDA). After performing PCR in the presence of another unlabeled primer, we separated the 362 bp PCR product with 2% low melting point agarose gel electrophoresis. The gel was then immersed into a Eu<sup>3+</sup> solution. During soaking, Eu<sup>3+</sup> diffuses into the gel and associates with BCPDA to form a fluorescent complex of long fluorescence lifetime. This complex can be quantified by scanning the gel with a time-resolved fluorometric reader. Because BCPDA and Eu<sup>3+</sup> are not fluorescent by themselves, background signals are very low. The detection limit was about 5 ng of DNA. We have also shown that the BCPDA-labeled product could be blotted and detected on the membrane by using an anti-BCPDA antibody. These two technologies may find applications other than in PCR, e.g. in fluorescence-based DNA sequencing and in solution hybridization.

PMID: 7679258, UI: 93158897

---

 the above report in   format documents on this page through Loansome Doc

---

NCBI PubMed Pub MEDLINE QUERY PubMed ?

Other Formats: Citation MEDLINE

Links: Related Articles

☐ Order this document*Clin Biochem* 1994 Oct;27(5):373-379

## Detection of Philadelphia chromosome using PCR and europium-labeled DNA probes.

Eskola JU, Hamalainen M, Nanto V, Rajamaki A, Dahlen P, Iitia A, Siitari H

Joint Clinical Biochemistry Laboratory, University of Turku, Turku University Hospital, Finland.

More than 95% of the patients with chronic myelogenous leukemia (CML) carry translocations between protooncogene *abl* of chromosome 9 and *bcr* gene of chromosome 22, resulting in the Philadelphia chromosome (Ph1). After allogeneic bone marrow transplantation (BMT) it is important to detect possible residual malignant cells in CML patients. A new sensitive hybridization method combined with polymerase chain reaction (PCR), based on the detection of the europium (Eu3+) label by time-resolved fluorescence, was applied for the detection of Ph1 chromosome. Total RNA from 10(6) peripheral blood leukocytes was isolated by the acid guanidinium thiocyanate-phenol-chloroform extraction. After cDNA synthesis by reverse transcriptase, the PCR amplification (30 cycles) was carried out. In the detection phase two oligonucleotide probes were used in the hybridization reaction, one biotinylated (*bcr* gene, exon 2) and one (*abl* gene) labeled with Eu3+. The hybrids were collected in a streptavidin-coated microtitration well and the bound Eu3+ was measured in a time-resolved fluorometer. To assess the sensitivity of the method, different numbers of CML cell line K562 cells were mixed with 10(5) apparently normal human leukocytes. Five K562 cells/10(5) leukocytes could be detected. Six patients with CML confirmed by clinical and cytogenetic criteria were studied. Three of the patients underwent an allogeneic BTM 6-18 months before the investigation and all of them were Ph1-negative. The other three patients who were nontransplanted were positive as expected.

PMID: 7867215, UI: 95171612

 the above report in   format documents on this page through Loansome Doc

NCBI PubMed Pub MEDLINE QUERY PubMed ?

Other Formats: Citation MEDLINE

Links: Related Articles

☐ Order this document*J Clin Microbiol* 1993 Jul;31(7):1886-1891

## Detection of adenovirus in clinical specimens by polymerase chain reaction and liquid-phase hybridization quantitated by time-resolved fluorometry.

**Hierholzer JC, Halonen PE, Dahlen PO, Bingham PG, McDonough MM**

Respiratory and Enteric Viruses Branch, Centers for Disease Control, Atlanta, Georgia 30333.

In addition to tests for the group-specific hexon antigen of adenoviruses, adenoviruses can be detected in clinical specimens by hybridization assays utilizing the widely shared base sequences of the region of the hexon gene that codes for the group-reactive determinants. We have developed a liquid-phase hybridization system with biotin- and europium-labeled probes which are reacted after DNA amplification of a 161-bp region of the hexon gene and which are quantitated by time-resolved (TR) fluorometry in streptavidin-coated microtiter wells. Polymerase chain reaction (PCR)-TR fluorometry is not a rapid test in the usual sense, but it is highly useful for specimens with inherent toxicity or with low virus yield, such as organ minces and specimens obtained late in the course of an illness. In a survey of 103 specimens tested by this method, including urine, stool, and tissue suspensions, the agreement with the hexon-specific TR fluoroimmunoassay antigen test for positive specimens was 100% and the sensitivity compared with that of virus culture was 91%. The PCR-TR fluorometry system was also shown to be advantageous as a quantitative measure of PCR products.

PMID: 8349768, UI: 93352755

 the above report in   format documents on this page through Loansome Doc



NCBI PubMed Pub MEDLINE QUERY PubMed ?

Other Formats: Citation MEDLINE

Links: Related Articles

☐ Order this document*Clin Chem* 1993 Feb;39(2):196-201

## Europium(III) trisbipyridine cryptate label for time-resolved fluorescence detection of polymerase chain reaction products fixed on a solid support.

Lopez E, Chypre C, Alpha B, Mathis G

CIS Bio International, Laboratoire des Produits pour l'Analyse Medicale, Bagnols sur Ceze, France.

We report a method for the nonradioactive detection of specific amplified target DNA, based on a new long-lived fluorescent label, europium(III) trisbipyridine cryptate [TBP-(Eu<sup>3+</sup>)]. After a polymerase chain reaction with a set of nested primers, microtiter affinity-collected hybrids were recognized by a TBP(Eu<sup>3+</sup>)-labeled antibody. The particular structure of the label containing the energy transfer units (bipyridine moieties) enables direct fluorescence measurement on the solid phase. Moreover, its stability and long lifetime contribute to the enhancement of the signal-to-noise ratio. The usefulness of this new label is illustrated by the detection of human papillomavirus type 16 DNA in clinical smears.

PMID: 8381731, UI: 93161526

 the above report in   format documents on this page through Loansome Doc

NCBI PubMed Pub MEDLINE QUERY PubMed ?

Other Formats: Citation MEDLINE

Links: Related Articles

☐ Order this document*Nucleic Acids Res* 1991 Jan 11;19(1):109-116

## Post-PCR sterilization: development and application to an HIV-1 diagnostic assay.

Isaacs ST, Tessman JW, Metchette KC, Hearst JE, Cimino GD

HRI Research, Inc., Berkeley, CA 94710.

We have developed an effective post-PCR sterilization process and have applied the procedure to a diagnostic assay for HIV-1. The method, which is based on isopsoralen photochemistry, satisfies both the inactivation and hybridization requirements of a practical sterilization process. The key feature of the technique is the use of isopsoralen compounds which form covalent photochemical adducts with DNA. These covalent adducts prevent subsequent extension of previously amplified sequences (amplicons) by Taq polymerase. Isopsoralens have minimal inhibitory effect on the PCR, are activated by long wavelength ultraviolet light, provide sufficient numbers of covalent adducts to impart effective sterilization, modify the amplified sequence such that it remains single-stranded, and have little effect on subsequent hybridization. The sterilization procedure can be applied to a closed system and is suitable for use with commonly used detection formats. The photochemical sterilization protocol we have devised is an effective and pragmatic method for eliminating the amplicon carryover problem associated with the PCR. While the work described here is limited to HIV-1, proper use of the technique will relieve the concern associated with carryover for a wide variety of amplicons, especially in the clinical setting.

PMID: 2011490, UI: 91187605

 the above report in   format documents on this page through Loansome Doc

NCBI PubMed Pub? QUERY PubMed ?

Other Formats: Citation MEDLINE

Links: Related Articles

☐ Order this document*Nucleic Acids Res* 1991 Jan 11;19(1):99-107

## Post-PCR sterilization: a method to control carryover contamination for the polymerase chain reaction.

Cimino GD, Metchette KC, Tessman JW, Hearst JE, Isaacs ST

HRI Research, Inc., Berkeley, CA 94710.

We describe a photochemical procedure for the sterilization of polynucleotides that are created by the Polymerase Chain Reaction (PCR). The procedure is based upon the blockage of Taq DNA polymerase when it encounters a photochemically modified base in a polynucleotide strand. We have discovered reagents that can be added to a PCR reaction mixture prior to amplification and tolerate the thermal cycles of PCR, are photoactivated after amplification, and damage a PCR strand in a manner that, should the damaged strand be carried over into a new reaction vessel, prevent it from functioning as a template for the PCR. These reagents, which are isopsoralen derivatives that form cyclobutane adducts with pyrimidine bases, are shown to stop Taq polymerase under conditions appropriate for the PCR process. We show that effective sterilization of PCR products requires the use of these reagents at concentrations that are tailored to the length and sequence of the PCR product and the level of amplification of the PCR protocol.

PMID: 2011516, UI: 91187655

---

 the above report in   format documents on this page through Loansome Doc

---

Other Formats: [Citation](#) [MEDLINE](#)Links: [Related Articles](#)☐ Order this document*J Clin Microbiol* 1993 Sep;31(9):2356-2360

## Preventing false positives: quantitative evaluation of three protocols for inactivation of polymerase chain reaction amplification products.

Rys PN, Persing DH

Division of Clinical Microbiology, Mayo Clinic, Rochester, Minnesota 55905.

False-positive results because of carryover contamination by previously amplified nucleic acids are currently the greatest impediment to routine implementation of nucleic acid amplification protocols. We evaluated three methods for inactivation of a 156-bp *Borrelia burgdorferi* polymerase chain reaction (PCR) product: (i) post-PCR cross-linking with isopsoralen (IP), (ii) pre-PCR treatment of a dU-containing PCR product with uracil N-glycosylase (UNG), and (iii) post-PCR alkaline hydrolysis (primer hydrolysis) of PCR products synthesized by using primers containing 3' ribose residues. The sensitivities of the PCR performed under the conditions of each protocol were comparable. Inactivation of amplified DNA was highly efficient for all three protocols; the IP and UNG protocols eliminated at least to  $3 \times 10(9)$  copies of the product. The primer hydrolysis protocol varied in efficiency depending on the number and position of the 3' ribose residues, but inactivation ranged from  $10(4)$  to  $10(9)$  copies. We conclude that with some modifications, all three systems are effective for eliminating amplified DNA products. Routine implementation of at least one method should help to avoid false-positive results because of carryover contamination.

PMID: 8408555, UI: 94013388

---

 the above report in [Macintosh](#)  format documents on this page through Loansome Doc

---

NCBI PubMed PubMed ?

Other Formats: Citation MEDLINE

Links: Related Articles

☐ Order this document*J Clin Microbiol* 1996 Dec;34(12):3115-3119

## Use of AmpliWax to optimize amplicon sterilization by isopsoralen.

**De la Viuda M, Fille M, Ruiz J, Aslanzadeh J**

Department of Laboratory Medicine, University of Connecticut School of Medicine, Farmington 06030, USA.

The photochemical inactivation of amplicons by isopsoralen (IP-10) has been suggested as a possible means to prevent PCR carryover contamination. To evaluate the technique, serial dilutions of amplicons (10(11) to 10(3)) from the *Borrelia burgdorferi* OSP A gene were amplified in the presence of 0, 25, 50, and 100 micrograms of IP-10 per ml for 45 cycles. The PCR products were exposed to UV light for 15 min to activate IP-10 and sterilize the amplicons. One microliter of each sterilized sample was reamplified for an additional 45 cycles. The PCR products were then resolved in an agarose gel, blotted onto a nylon membrane, and probed with an alkaline phosphatase-conjugated chemiluminescent probe. Although IP-10 at concentrations of 50 and 100 micrograms/ml effectively sterilized up to 10(11) amplicons, the compound was inhibitory to PCR. IP-10 at a concentration of 25 micrograms/ml had slight inhibitory effect on PCR and did not completely sterilized all of the amplicons. Therefore, in subsequent experiments AmpliWax was substituted for mineral oil, and PCR was performed on 10(9) to 10(3) amplicons as described above. Following the amplification, the PCR tubes were cooled to solidify the AmpliWax and inoculated with various concentrations of IP-10. With this technique, PCR products produced from as many as 10(9) target amplicons were effectively sterilized with 200 micrograms of IP-10 per ml. Similarly, the addition of IP-10 (50 micrograms/ml) before and after PCR was evaluated for the detection of *B. burgdorferi* in 62 ticks from a region of Southern Connecticut where the organism is highly endemic. PCR performed in the presence of 50 micrograms of IP-10 per ml detected *B. burgdorferi*-specific DNA in 17 of 62 ticks (27%) following gel electrophoresis and in 34 of 62 ticks (55%) following Southern blot hybridization of the PCR products. In contrast, post-PCR addition of IP-10 detected borrelia-specific DNA in 31 of 62 ticks (50%) following gel electrophoresis and in 46 of 62 ticks (64%) following Southern blot hybridization. We conclude that the replacement of mineral oil with AmpliWax can be useful in eliminating the inhibitory effects of IP-10 and other sterilizing agents for post-PCR sterilization of amplicons.

PMID: 8940457, UI: 97095036

 the above report in   format documents on this page through Loansome Doc

National Library  
of Medicine

PubMed

[PubMed](#)[Nucleotide](#)[Protein](#)[Genome](#)[Structure](#)[PopSet](#)[Taxonomy](#)[OMIM](#)Search 

for

[Limits](#)[Preview/Index](#)[History](#)[Clipboard](#)[About Entrez](#)[Entrez PubMed](#)[Overview](#)[Help | FAQ](#)[New/Noteworthy](#)[PubMed Services](#)[Journal Browser](#)[MeSH Browser](#)[Single Citation Matcher](#)[Batch Citation Matcher](#)[Clinical Queries](#)[Cubby](#)[Related Resources](#)[Order Documents](#)[Grateful Med](#)[Consumer Health](#)[Clinical Alerts](#)[ClinicalTrials.gov](#)■ 1: *J Protozool* 1989 Jan-Feb;36(1):60S-61S[Related Articles, Books](#)

## Analysis of a surface antigen of *Pneumocystis carinii*.

**Linke MJ, Walzer PD**

VA Medical Center, Cincinnati, OH.

The 115 kd band in polyacrylamide gels is a major antigen of *Pneumocystis carinii*. Data obtained from treatment with enzymes, binding to lectins, and labelling the surface with biotin suggest that this moiety is a glycoprotein containing mannosyl/glucosyl and N-acetylglucosamine residues, and that it is located on the cell wall of the organism. Other rat and human *P. carinii* antigens also are glycoproteins but differ in specific protein or carbohydrate residues or location on the organism.

PMID: 2785200

[Privacy Policy](#)[Write to the Help Desk](#)[NCBI | NLM | NIH](#)[Department of Health & Human Services](#)[Freedom of Information Act | Disclaimer](#)